

## FastClix Lancing Device + Accu-chek Guide Meter



Set-up:  
Turn the cap to the number you want.  
1 is for sensitive skin  
5 is for thicker skin  
2-3-4 are inbetween



1. Pull cap straight off



2. Push drum in (each drum has 6 lancets)




3. Make sure drum is all the way in



4. Replace cap until it clicks back on. "6" new lancets can be used



5. Remove one test strip

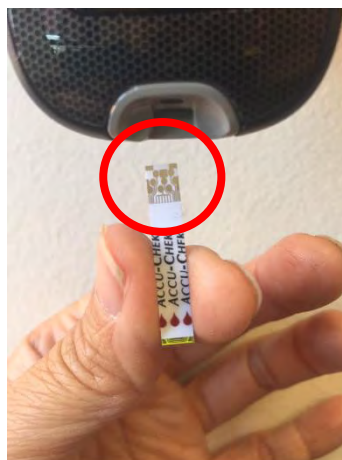

6. Make sure gold side goes into meter



7. Insert test strip. No need to push any buttons.



7. Push strip all the way in. Meter will turn on automatically




8. Touch tip of lancing device to side of finger



9. Push button (like a ballpoint pen) to poke your finger



10. Squeeze finger if you need a bigger drop of blood

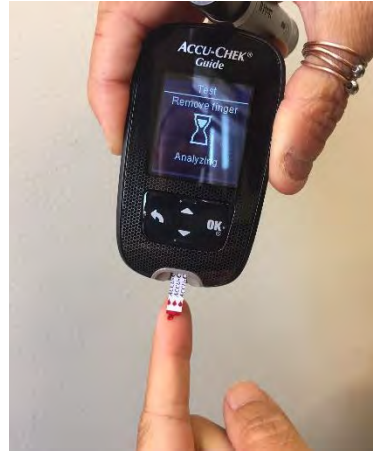

11. Apply blood anywhere on yellow end of test strip



12. After getting your result, pull test strip out to turn meter off



13. Move to next lancet by switching the white tab over and back



14. Lancet counter will go down one number. Change the whole drum when you get to "1"

Replace drum after 6 lancets are used up:



Remove cap



Pull out lancet drum



Throw used drum away into sharps box

### Other important reminders

- There are 6 lancets in each lancet drum. Make sure you use them all before taking out the drum
- Before checking your sugar, wash your hands with soap and water
- Don't share your lancing device with other people
- A battery symbol at the top of the meter screen or "E-9" error message means it's time to change the battery
- Check the expiration on the bottle of test strips. Don't use expired test strips

